GRADIA® DIRECT

LIGHT-CURED COMPOSITE RESTORATIVE

For use only by a dental professional in the recommended indications.

RECOMMENDED INDICATIONS

1. Direct restorative for Classes III, IV, V cavities.
2. Direct restorative for wedge-shaped defects and root surface cavities.
3. Direct restorative for veneers and diastema closure.
4. Direct restorative DIRECT POSTERIOR (radioopaque)
5. Direct restorative for Class I and II cavities

CONTRAINDICATIONS

1. Pulp capping.
2. In rare cases the product may cause sensitivity in some people. If any such symptoms are experienced, discontinue the use of the product and refer to a physician.

DIRECTIONS FOR USE


1. Shade Selection
Clean the tooth with pumice and water. Shade selection should be made prior to isolation. Select the appropriate GRADIA DIRECT shades by referring to the GRADIA DIRECT shade guide.
2. Cavity Preparation
Prepare cavity using standard techniques.
- Note:
For pulp capping, use calcium hydroxide.
3. Bonding Treatment
For bonding GRADIA DIRECT to enamel and / or dentin, use a light-cured bonding system such as GC Fuji BOND LC, GC UniFil® Bond or GC G-BOND™ (Fig. 1). Follow manufacturer's instructions.

4. Placement of GRADIA DIRECT
1) Dispensing from a Unitip
Insert the GRADIA DIRECT Unitip into the Unitip APPLIER or equivalent. Remove the cap and extrude material directly into the prepared cavity. Use steady pressure (Fig. 2). Maintain pressure on the applicator handle while removing the Unitip APPLIER and Unitip from the mouth. This prevents the Unitip from coming loose from the applicator.
- 2) Dispensing from a syringe
Remove syringe cap and dispense material onto a mixing pad. Place the material into the cavity using a suitable placement instrument. After dispensing, screw syringe plunger anticlockwise by a half to full turn to release residual pressure inside the syringe. Replace cap immediately after use.
- Note:
1. Basically, material can be applied in a single layer to achieve aesthetic restorations using Standard shades. For details, refer to the Clinical Hints.
2. Material may be hard to extrude immediately after removing from cold storage. Prior to use, leave to stand for a few minutes at normal room temperature.
3. After dispensing, avoid too long an exposure to ambient light. Ambient light can shorten the manipulation time.

5. Clinical Hints
a. Anterior cavities
In the case of small cavities
Restore using a one shade layering technique. In most cases the use of one Standard shade alone will be sufficient. In cases where a higher degree of translucency is needed, one of the Outside special shades can be selected. See also Examples of Clinical Applications.
- b. In the case of large cavities
In most cases a multi shade layering technique will give the best aesthetic results. To block out shine throughs from the oral cavity or to mask discoloured dentin, select an appropriate Inside special shade and continue to build up with a Standard shade.
To make a restoration more life like e.g. to copy age related changes in appearance, the final layer should consist of an Outside special shade. See also Examples of Clinical Applications and / or consult the Shade Combination Chart.
6. Posterior cavities
a. In the case of small cavities
Restore using a one shade layering technique. In most cases the use of one Standard shade alone will be sufficient. In cases where a higher translucency is needed, one of the Outside special shades should be selected. See also Examples of Clinical Applications.
- b. In the case of deep cavities
Place a flowable composite such as GRADIA DIRECT Flo or GRADIA DIRECT LoFlo™ on the cavity floor. Then place a Standard shade. For optimal aesthetics use an Outside special shade as the final composite layer. See also Examples of clinical applications.

*GC Fuji LINING® PASTE PAK, GC Fuji LINING® LC or GC Fuji IX GP can also be used as a liner or base material. Follow the respective manufacturer's instructions for use.

Examples of Clinical Applications (Clinical Hint No. 1, 2)

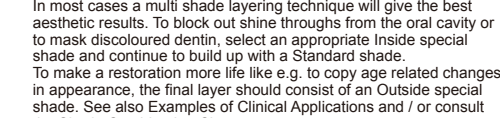Shade combination chart for multiple layers in large cavities of anterior teeth (1-b)

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For details of shades, refer to the following section of SHADES.

5. Contouring before Light Curing
Contour using standard techniques.
6. Light Curing
Light cure GRADIA DIRECT using a light curing unit (Fig. 3).
Refer to the following chart for Irradiation Time and Effective Depth of Cure.

GRADIA DIRECT ANTERIOR: Irradiation Time and Effective Depth of Cure

Irradiation time	Plasma arc (2000mW/cm²)	3 sec.	6 sec.
	GC G-Light (1200mW/cm²)	10 sec.	20 sec.
	Halgén / LED (700mW/cm²)	20 sec.	40 sec.
CT, NT, WT, GT, CVT		3.0 mm	3.5 mm
A1, A2, B1, B2, XBW, BW, DT		2.5 mm	3.0 mm
A3, B3		2.0 mm	3.0 mm
A3.5, A4, C3, A02, A03, A04, CV, CVD		1.5 mm	2.5 mm

GRADIA DIRECT POSTERIOR: Irradiation Time and Effective Depth of Cure

Control using standard techniques.

6. Light Curing

Light cure GRADIA DIRECT using a light curing unit (Fig. 3).

Refer to the following chart for Irradiation Time and Effective Depth of Cure.

GRADIA DIRECT ANTERIOR: Irradiation Time and Effective Depth of Cure

Irradiation time	Plasma arc (2000mW/cm²)	3 sec.	6 sec.
	GC G-Light (1200mW/cm²)	10 sec.	20 sec.
	Halgén / LED (700mW/cm²)	20 sec.	40 sec.
P-A1, P-A2, P-WT, P-NT		2.5 mm	3.0 mm
P-A3, P-A3.5		2.0 mm	3.0 mm

Fig. 3

- Note:
1. Material should be placed and light cured in layers. For maximum layer thickness, please consult the tables.
2. Lower light intensity may cause insufficient curing or discoloration of the material.
7. Finishing and Polishing
Finish and polish using diamond burs, polishing points and discs. To obtain a high gloss, polishing pastes can be used.

Cleaning and Sterilizing of the Unitip APPLIER

1. Clean the applicator with gauze soaked in alcohol.
2. The applicator may be steam autoclaved at 121°-135°C (250-275°F) for 10 to 30 minutes or disinfected using alcohol.
- Note:
1) Do not attempt to sterilize in a unit at a temperature that exceeds 135°C (275°F) during the drying step of the autoclave procedure. Otherwise, the applicer may be damaged.
- 2) The use of various chemical disinfectants / sterilants may cause damage to the applicator and therefore is not recommended.

Cleaning of the Shade Guide

The shade guide can be cleaned with water and mild soap. It should not be sterilized. The use of various chemical disinfectants / sterilants may cause damage to the shade guide.

SHADES

- 1.22 shades for anteriors
Standard shade : XBW (Extra Bleaching White), BW (Bleaching White), A1, A2, A3, A3.5, A4, B1, B2, B3, C3, CV (Cervical), CVD (Cervical Dark), CVD (Cervical Dark)
- Outside special shades : WT (White Translucent), DT (Dark Translucent), CT (Clear Translucent), NT (Natural Translucent), GT (Gray Translucent), CVT (Cervical Translucent)
- 2.6 shades for posteriors
Standard shade : P-A1, P-A2, P-A3, P-A3.5, Outside special shades : P-WT (White Translucent), P-NT (Natural Translucent)
- Note:
A, B, C, AO shades are based on Vita®* Shade.

STORAGE

Recommended for optimal performance, store in a cool and dark place (4° - 25°C / 39.2° - 77.0°F) away from high temperatures or direct sunlight. (Shelf life)

GRADIA DIRECT ANTERIOR - 3 years from date of manufacture

GRADIA DIRECT POSTERIOR - 3 years from date of manufacture

PACKAGES

1. Unitips
1. INTRODUCTION KIT (35 Unitips in 7 shades)
- 35 tips (5 tips each in 7 shades) (0.16ml per tip)
(6 anterior shades : A2, A3, A3.5, A03, CV, CT)
(1 posterior shade : P-A2)
- Shade guide
2. REFILL
- 1 syringe (in 28 shades) (2.7ml per syringe)
- Note:
Weight per Unitip : 0.24g for anterior tip, 0.28g for posterior tip
3. Option
a. Unitip APPLIER b. Shade guide c. Mixing pad (No.14B)

2. Syringes
1. INTRODUCTION KIT (7 syringes in 7 colors)
(6 anterior shades : A2, A3, A3.5, A03, CV, NT, DT)
(1 posterior shade : P-A2)
- Shade guide
2. REFILL
- 1 syringe (in 28 shades) (2.7ml per syringe)
- Note:
Weight per syringe : 4.0g for anterior syringe, 4.7g for posterior syringe
3. Option
a. Shade guide b. Mixing pad (No.14B)

CAUTION

1. In case of contact with oral tissue or skin, remove immediately with cotton or a sponge soaked in alcohol. Flush with water.
2. In case of contact with eyes, flush immediately with water and seek medical attention.
3. Take care to avoid ingestion of the material.
4. Wear plastic or rubber gloves during operation to avoid direct contact with air inhibited resin layers in order to prevent possible sensitivity.
5. For infection control reasons, Unitips are for single use only.
6. Wear protective eye glasses during light curing.
7. When polishing the polymerised material, use a dust collector and wear a dust mask to avoid inhalation of cutting dust.
8. Do not mix with other similar products.
9. Avoid getting material on clothing.
10. In case of contact with unintended areas of tooth or prosthetic appliances, remove with instrument, sponge or cotton pellet before light curing.
11. Do not use GRADIA DIRECT in combination with eugenol containing materials as eugenol may hinder GRADIA DIRECT from setting.

Last revised : 01/2015

MANUFACTURED BY
GC DENTAL PRODUCTS CORP.
2-289 Tomimachi-cho, Kanagawa, Japan 480-0840, Japan
RESPONSIBLE MANUFACTURER IN CANADA:
GC AMERICA INC.
3737 West 127th Street, Mississauga, ON, M5B 0B3 U.S.A.

DISTRIBUTED BY
GC CORPORATION
761-1 Heumachi-cho, Bishu-shi, Tokyo 174-8585, Japan
EU - GC EUROPE N.V.
R-3001 Leuven, Belgium - T: +32 18 10 10 00

GC AMERICA INC.
3737 West 127th Street, Mississauga, ON, M5B 0B3 U.S.A.
TEL: +1 905 587-5900

GC ASIA DENTAL PTE. LTD.
11 Tanjong Pagar Road, Singapore 068279
TEL: +65 6546 7588

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