

DESENSITISING VARNISH WITH FLUORIDES



HIGHLY EFFECTIVE

Hypersensitivity

Sensitive teeth are a widespread problem for a large number of people. Usually, the condition causes only transient pain which patients do not always report to their dentists. Hypersensitivity is mainly caused by exposed dentinal tubules. Therefore, the aim of treatment is to seal these tubules.

Indications

Profluorid Varnish is a fluoride-containing varnish for treating hypersensitivity. Dentinal tubules exposed on the necks of teeth or after cavity preparation are effectively sealed.


Components and effect of VOCO Profluorid® Varnish

Profluorid Varnish is a colophony-based varnish containing 5 % sodium fluoride (22,600 ppm fluoride). The fluoride ion, together with the calcium ions accumulated in the tubules, causes a precipitation of calcium fluoride, which effectively seals the tubules. According to a user survey, over 85 % of patients were completely free of pain after only one treatment with Profluorid Varnish.



As well as rapidly sealing the dentinal tubules, Profluorid Varnish also causes calcium fluoride deposits to form on the tooth surface. These deposits protect the tooth from acid attacks, promote remineralisation and contribute to the long-term formation of fluoroapatite.

Xylitol, which has a proven cariostatic effect in addition to its taste-enhancing properties, is also added to the varnish.

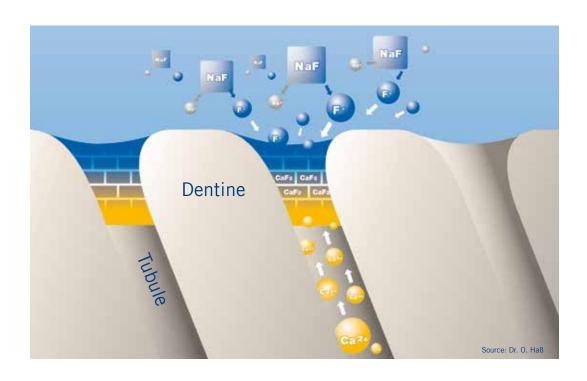


SIMPLE AND AESTHETIC

Easy to apply

Thanks to its colophony matrix, Profluorid Varnish adheres well even to wet surfaces. It is not necessary to dry the tooth surfaces. This improves cooperation of patients with low levels of compliance (e.g. children).

Aesthetics and taste

Profluorid Varnish is a tooth-shaded varnish which does not impair aesthetics. Profluorid Varnish is available in four pleasant flavours: melon and cherry – aimed especially at children and younger patients, mint – a classic, particularly for adults, and caramel – to provide a slightly more exotic option. These four refined flavours will enable you to provide an even better service by tailoring it to the individual patient.

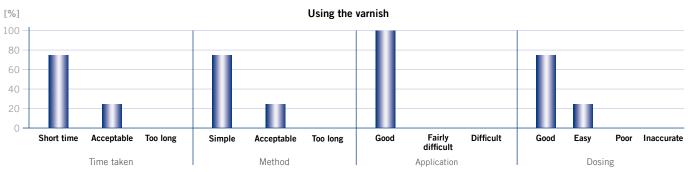

Survey of testing dentists

Opinion of dentists

Dentists using the product all confirmed its good handling and aesthetic properties. A vast majority also reported positive differences between Profluorid Varnish and competing products. Profluorid Varnish was equally well accepted by the patients.

After the treatment, nearly all of the patients reported that they were either symptom-free or had a significantly reduced level of sensitivity.

Profluorid Varnish: much more than just good taste!

Simple and aesthetic desensitising with VOCO Profluorid® Varnish







The SingleDose advantages
Profluorid Varnish SingleDose
comes in a blister pack containing
a brush and the varnish. This
enables the product to be applied
under perfectly hygienic conditions.

Perfect application with VOCO Profluorid® Varnish







The cartridge advantages The tried and tested cartridge enables you to apply the material with pinpoint accuracy, and to fluoridate a complete dental arch in one go.

YOUR CHOICE

Indications

Treatment of hypersensitive teeth

Sealing of the dentinal tubules on sensitive root surfaces or following cavity preparation



Advantages

- Easy to apply on wet surfaces
- · Aesthetic, white-transparent tooth colour
- Fast desensitisation and fluoride release
 (5 % NaF ≜ 22,600 ppm fluoride)
- Available in four flavours: melon, mint, cherry and caramel

Advantages of the SingleDose

- Product and applicator in hygienic disposable packing
- No additional devices required
- No spilling in any position

Advantages of the cartridge

- Suitable for multiple use with disposable cannulae
- Perfectly precise application
- Extremely easy to apply in proximal areas
- Fluoridates a complete dental arch in one go

Presentation

50 ml REF 1272	5 × 1.70 ml REF 2235
REF 1272	REF 2235
_	_
_	_
_	_
_	_
	-

VOCO GmbH
Anton-Flettner-Straße 1-3
27472 Cuxhaven
Germany

Tel.: +49 (0) 4721-719-0 Fax: +49 (0) 4721-719-140

info@voco.com www.voco.com

Available from:		

